

Tissue Homogenizer

WISBIOMED

Tel: 650-242-8368
Fax: 650-242-8244

Homogenize 24 samples in two minutes



An efficient tissue homogenizer based on a new technology is being offered at a low price.


Works on plants, animal tissue, organ samples, and cells.

From tree leaves to cultured cells, all can be homogenized in a few minutes with the bead based homogenizer

Easy to operate:

1. Place tissue pieces in a tube and add beads and buffer
2. Place tube in the platform
3. Close lid, set time and speed, and push start
4. Take the tubes out

Process 24 samples with small volumes. Or 9 samples with large volumes.

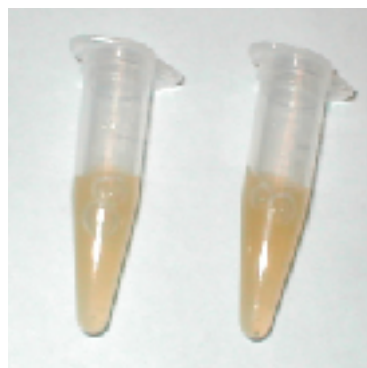

Homogenized solution looks hazy & uniform
No chunks.

Available beads:

1. Bulk package. different sizes and densities
2. Pre-loaded tubes with ZrO beads

Visit www.wisbiomed.com
for more information.

